


Rocket® BLUE™

Multipurpose Needle Drain

Rocket® BLUE™ Multipurpose Needle Drain offer a purpose designed, cost-effective solution for routine fluid drainage procedures in the department or ward settings.

Indications:

For the drainage of body fluids including ascites, urine and pleural fluid. For insertion into, but not restricted to, the abdominal cavity, chest and urinary bladder using a 'direct stab' technique combined with aspiration and/or ultrasound guidance.

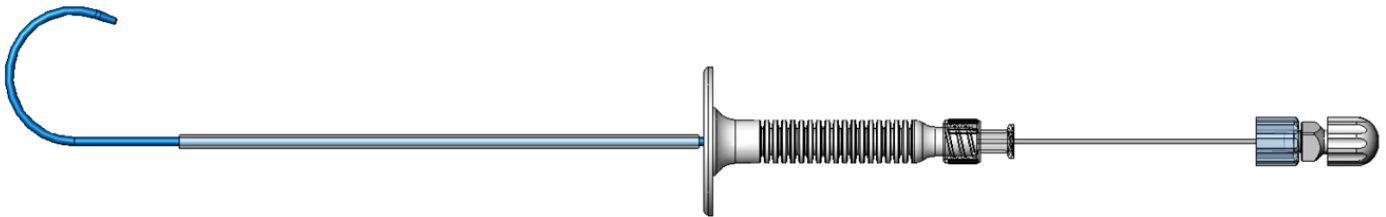

- **CE MARKED, PURPOSE DESIGNED DEVICE**, designed specifically for paracentesis and drainage of abdominal ascites, urinary bladder and other drainage applications which avoids the risks of inappropriate use of non specific devices such as generic suprapubic drains.
- **SOFT, KINK FREE 6FG** catheter material, improves patient comfort and prevents unnecessary catheter blockage caused by thin walled materials.
- **1CM CATHETER MARKINGS** to aid safer, controlled insertion.
- **SOFT INTEGRAL STITCH PLATE** combined with flexible tube support, provides simple, secure fixation after insertion.
- **REMOVABLE NEEDLE CAP** permits aspiration to aid correct placement during insertion.
- **DRAINAGE BAG CONNECTION SET** included to make insertion and connection to a wide range of drainage options quick and secure.

Description	Order Code	NHS Code
<p>Rocket® BLUE™ Multipurpose Needle Drain Set 6FG x 20cm x 6 holes with integral stitch plate and insertion needle set. 30cm drainage bag connection set with 2 way stopcock, No.11 scalpel and Instructions for Use. Sterile, for single patient use, supplied in multiples of 5 sets.</p>	R58800-06-BD	FSW486

